

A close-up photograph of a woman's face, focusing on her mouth and nose. She has dark, glossy lips and a slight smile. Her hands are visible in the foreground, with long, manicured nails. The background is a soft, out-of-focus light color.

Platinum 3 in 1

LASER FOR SKIN REJUVENATION,
HAIR REMOVAL AND TATTOOS

Platinum 3 in 1

LASER



WHY PLATINUM ?

■ DIODE + ND YAG LASERS + Q SWITCH

- 1 Expand your treatment offerings for clients
- 2 Multi Platform Laser system
- 3 Efficiently reduce costs and conserve space, featuring a built-in cooling system for superior results and client comfort

■ OUR SUPPORT

- 4 Comprehensive Training Provided
- 5 Rent-to-Buy options available
- 6 Two-Year Limited Warranty
- 7 Product Launch and Marketing Support

■ MACHINE QUALITY

- 8 Italian Water pump for Better Cooling
- 9 Laser Module with German Dilas Bars
- 10 TDK 6 Ways Japanese Power Supply





PLATINUM 3 IN 1

TRIPLE WAVELENGTH DIODE LASER

Diode laser hair removal comes with a 755 + 808 + 1064 triple wavelength laser stack with semiconductors technology.

They are assembled to produce light for permanent hair removal, hair reduction, depilation, and epilation treatments



3 Wavelengths in 1 handpiece for Superior Results

Alex 755nm

The Alexandrite wavelength improves the absorption of powerful energy by the melanin chromophore, making it ideal for the most extensive range of hair colors and types, especially light and thin. With increased superficial penetration, the 755nm wavelength targets the Bulge of the hair follicle and is especially effective for superficially embedded hair.



Speed 808nm


The 808nm is a classic hair removal wavelength, offering deep penetration of the hair follicle with high average power, and a rapid repetition rate, and a large 2 cm²—spot size for time-efficient treatments. The 808nm has a moderate melanin absorption level, making it safe for darker skin types. Its deep penetration capabilities target the Bulge and Bulb of the hair follicle, while moderate tissue depth penetration makes it ideal for treating the arms, legs, cheeks and chin.



YAG 1064nm

The YAG 1064 wavelength is characterized by lower melanin absorption, making it a perfect solution for darker skin types. The 1064nm also offers the deepest penetration of the hair follicle, so it targets the Bulb and Papilla while deeply treating embedded hair in areas such as the scalp, armpits, and pubic areas. Higher water absorption generates a higher temperature, incorporating the 1064nm wavelength increases the thermal profile of the overall treatment, leading to highly effective hair removal and improved results

Most Advanced Hair Removal Technology

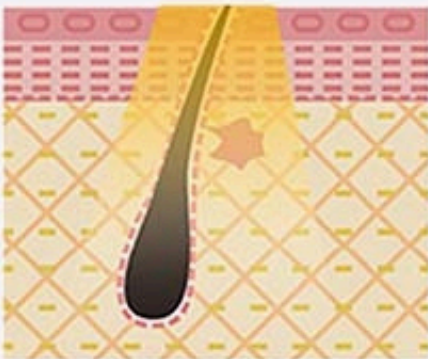

Emit a laser beam




The laser beam is absorbed by hair follicle melanin



Different color hair has different ability to absorb light wave



After absorption, the laser energy will be converted into heat energy



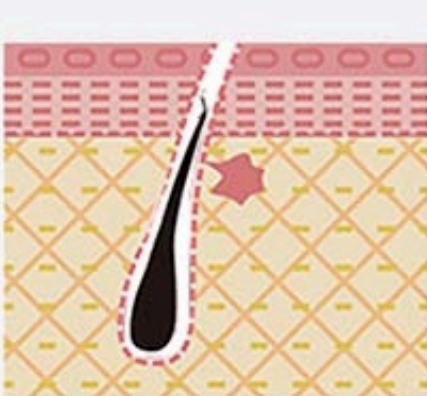


Hair follicles and surrounding stem cells lose their activity and fall off



So reach hair removal effect

Platinum 3 in 1 can treat any Body Area

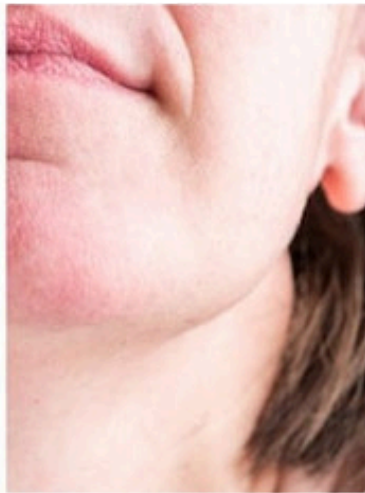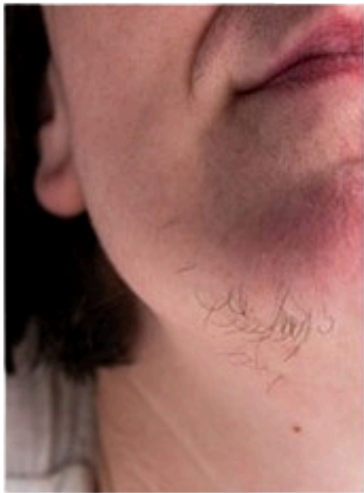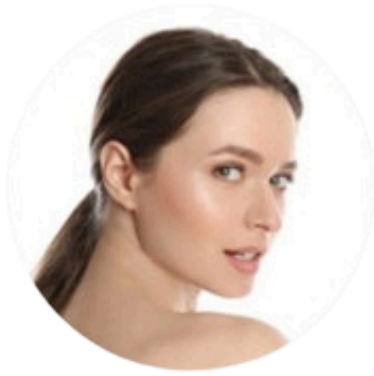

Hair Removal - Protocols

After thorough consultation begin the treatment following these steps.

- Clean the treatment area with a non-alcoholic antiseptic cleanser.
- Remove hair by shaving from the treatment area (advise the client to come shaved for the treatment).
- Clean the residue of the shaved hair which are present on the skin.
- Take a picture of the treatment area.
- Apply conductive gel to the treatment area.
- Begin laser hair removal treatment by selecting accurate parameters according to the skin type.
- Glide the handpiece on the treatment area in vertical, horizontal, and circular motions, giving 5 to 7 passes
- The tip of the handpiece must be perpendicular to the skin, and a degree of pressure must be applied to the skin. You should move the tip of the handpiece at a constant speed.
- Faster speed is used for the larger areas. One area should be 10 x 5 cm.
- The treatment area should be stretched so that there is no creasing, to avoid any epidermal damage.

Steps of the Treatment



Use a new shaving razor



Take a photo (with Consent)



Shave the Treatment Area



Remove the Residual Hair



Apply a Light Layer of Conducting Gel



**Disinfect the Laser Tip with
non alcoholic wipe**



Perform Laser Treatment



**Clean Treatment Area and
Apply Repair Gel**



Sanitize Laser Head / Tip

**Wear Specific Eye Protection
for everyone present in the
Treatment Room**





Make sure laser handpiece is fully in contact with the Skin and is not tilted or lifted




- The parameters are reference only, it is recommended to start with the low energy and adjust the parameters appropriately according to the customer's feedback during the treatment.
- The specific number of Light Spots / Laser Shots is for reference only it will be different according to the physical condition of different customers and at the endpoint response of the treatment shall prevail.

After the Treatment



It is normal to feel a slight burning / Tingling sensation after the treatment.

Mild and small papules will appear on the skin on hair follicles.
Please refer to the pictures above.

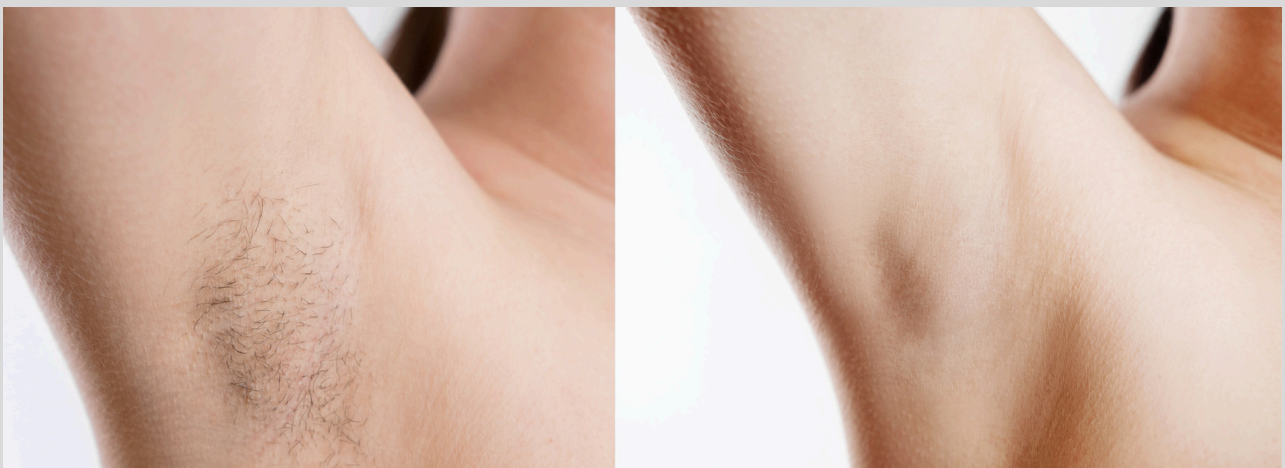

Treatment Areas

For smaller treatment areas e.g. lips, armpits & hairline give three to nine shots at the same point and cover the whole area.

For Larger areas e.g. limbs, back, abdomen, and legs, other sliding operations can be performed in multiple directions as described earlier, allowing the light/laser energy to better penetrate to the hair follicle.

For large areas, you have the option to use the foot pedal if preferred

						
	TYPE I	TYPE II	TYPE III	TYPE IV	TYPE V	TYPE VI
BEFORE SUN	Ivory	Fair or pale	Fair to beige, with golden undertones	Olive or light brown	Dark brown	Deeply pigmented dark brown to darkest brown
AFTER SUN	Always freckles, always burns/peels, never tans	Usually freckles, often burns/peels, rarely tans	Might freckle, burns on occasion, sometimes tans	Doesn't really freckle, rarely burns, often tans	Rarely freckles, almost never burns, always tans	Never freckles, never burns, always tans



I



II



III



IV



V



VI

Skin Type I



Hair Colour

Energy
 J/cm^2

Frequency
Hz

Coarse Black
Hair

40-65

8-10

Medium Dark
Brown Hair

45-65

9-12

Light Brown
Hair

50-70

10-16

Skin Type II



Hair Colour

Energy
 J/cm^2

Frequency
Hz

Coarse Black
Hair

40-60

8 - 12

Medium Dark
Brown Hair

45-65

9-12

Light Brown
Hair

45-70

10-16

Skin Type III



Hair Colour

Energy
 J/cm^2

Frequency
Hz

Coarse Black
Hair

40-60

8 -10

Medium Dark
Brown Hair

45-65

8-12

Light Brown
Hair

50-70

8-14

Skin Type IV



Hair Colour

Energy
 J/cm^2

Frequency
Hz

Coarse Black
Hair

40-60

8 - 10

Medium Dark
Brown Hair

45-65

8 - 12

Light Brown
Hair

50-70

8 - 14

Skin Type V



Hair Colour

Energy
 J/cm^2

Frequency
Hz

Coarse Black
Hair

35-45

7 - 9

Medium Dark
Brown Hair

45-49

7 - 9

Light Brown
Hair

45-55

7 - 9

Skin Type VI



Hair Colour

Energy
 J/cm^2

Frequency
Hz

Coarse Black
Hair

35-45

6 - 9

Medium Dark
Brown Hair

45-49

6 - 9

Light Brown
Hair

50-55

6 - 9

Parameters for Small Area

Energy from 8-12 J/cm² , 5-9 Hz



Number of Shots by Treatment Area

Face

350 - 650 Shots / Session

Each Armpit

300 - 400 Shots / Session

Abdomen / Back

1200- 2500 Shots / Session

Bikini

500 - 650 Shots / Session

Full Legs


1500- 4000 Shots / Session

Full Arms

1500- 3000 Shots / Session

OFFER YOUR CLIENTS

Carbon Laser Peel

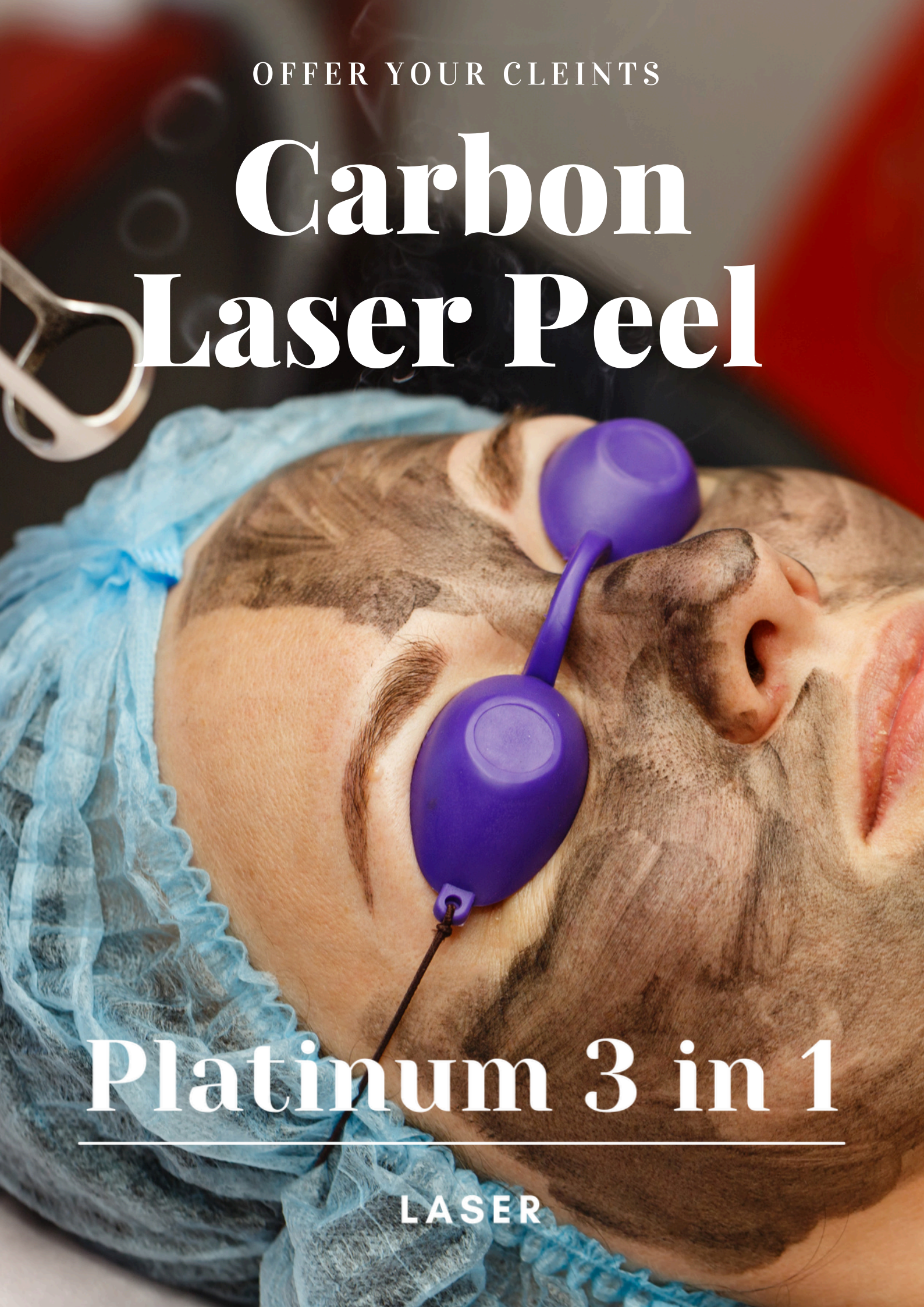

Platinum 3 in 1

LASER



Bring the skin treatment of
celebrities to your clients
with Platinum 3 in 1.

Carbon Peel is a non-invasive, completely painless laser-based skin treatment that quickly and effectively treats a variety of skin imperfections, with no side effects or patient downtime and with immediate results.

Using formulated unique carbon liquid applied to the face, in conjunction with the powerful and versatile laser, the treatment exfoliates the skin, unclogs and tightens pores, reduces oil production, and smooths rough and uneven skin texture.

With a unique 2-stage approach, Platinum Laser maximizes the effectiveness of the carbon peel with two distinct modes of laser delivery, while a special applicator allows aestheticians to treat the area from a distance, for better visibility, quick and a clean procedure.

The treatment is suitable for all skin types and is particularly beneficial for clients with oily skin, blackheads, enlarged pores, dull skin, rough/uneven skin texture, acne, fine lines, and wrinkles.

Two-stage approach

During the first stage, the laser gently heats the water molecules in the cream, allowing the carbon to bond well to the skin and to the contents within the pores.

At the second stage, the short pulse of the laser breaks up the carbon into fine particles. The laser wavelength is highly absorbed by the particles, thereby destroying the carbon and removing the dead skin cells, oils and impurities that were bound to it.

Procedure for the Carbon Laser Peel

- We recommend that the client stops using retinol for around one week before the laser carbon peel. They should also use sunscreen daily during this time.
- The laser carbon peel is a multi-part process that takes around 45 minutes from start to finish. For this reason, it is sometimes referred to as a lunchtime peel.
- Clean the face / Treatment area thoroughly.
- Spread a layer of liquid carbon on the treatment area. The carbon absorbs oil, dirt, and dead skin cells from deep within the pores.
- Wait for about 10 minutes for the liquid carbon to dry.
- Use protective eyewear for everyone present in the treatment room
- Pass the laser over the treatment area according to the skin in a low-energy setting. This heats the carbon, binding it to the epidermis.
- If the client has very sensitive skin or dark patches like melasma, your doctor may skip the heating step. The laser also heats subdermal tissue, boosting collagen and plumping up fine lines.
- On a higher-energy setting, use the laser to fragment the liquid carbon into tiny particles.
- Apply a hydrating moisturizer and sunscreen to the skin after treatment.

If the client has sensitive skin, they may experience a slight pinking or reddening of the skin. This usually lasts for one hour or less.

Parameters

TREATMENT NUMBER	HZ	ENERGY
1	1	600
2	2	700
3	3	700
4	3	700
5	3	800
6	4	800



Skin Rejuvenation

PLATINUM
SkinTech

Advanced Laser Technology for Clearer, Smoother And Healthier Skin

Platinum 3 in 1 Skin Rejuvenation treatments are designed to improve overall skin quality by stimulating collagen production, refining pores, reducing pigmentation, and enhancing skin tone – all with minimal downtime.

SAFE FOR ALL SKIN TYPES

The Nd:YAG laser technology used in Platinum 3 in 1 is known for its low melanin absorption, making it safe and effective for all skin types, including darker skin tones. This allows clinics to confidently treat a diverse client base while maintaining excellent safety and comfort standards.

Ideal for clients with:



Dull skin



Acne-prone skin




Oily skin



Fine lines



Pigmentation

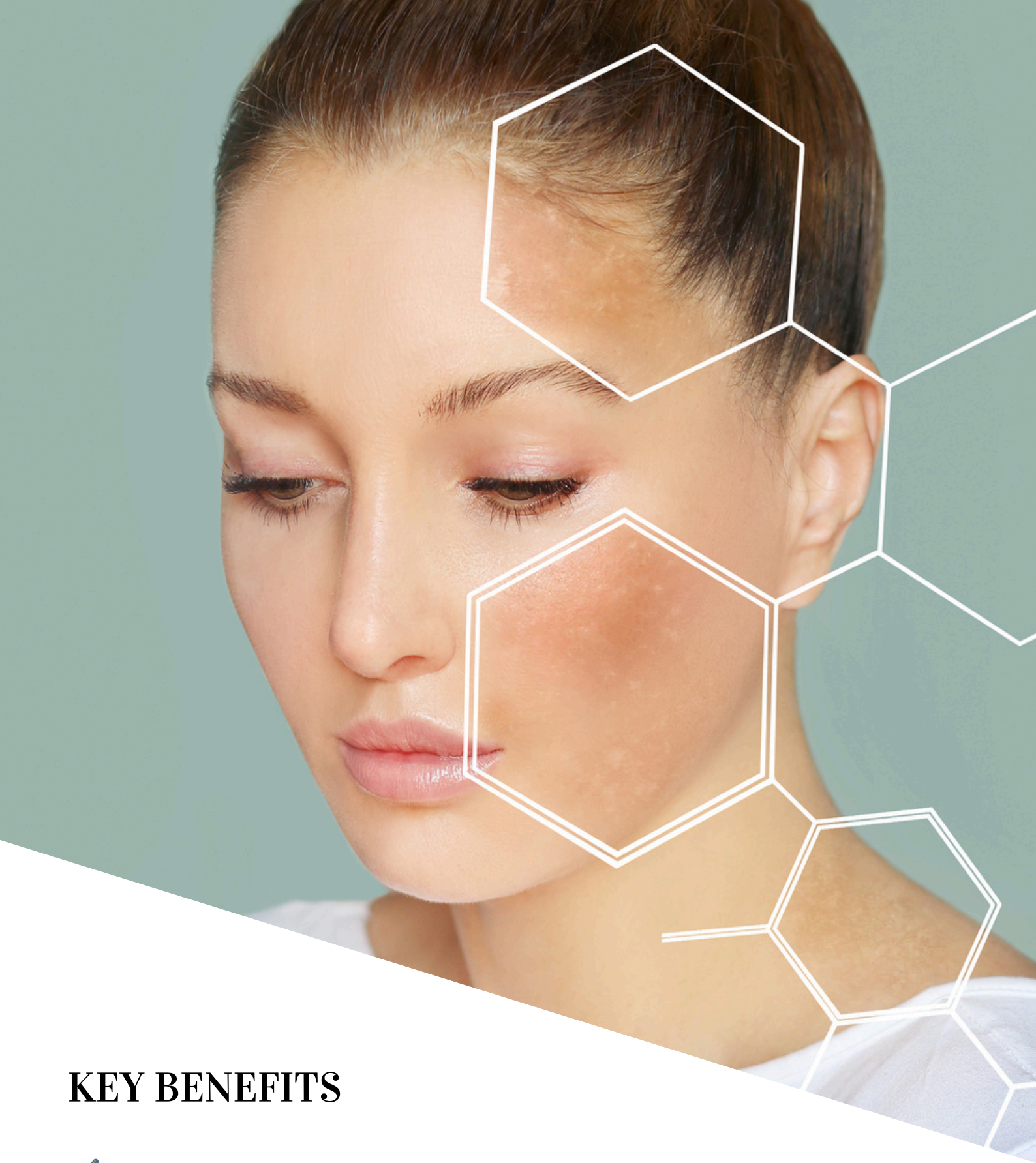

KEY BENEFITS

- ✓ Improves skin texture and tone
- ✓ Reduces fine lines and enlarged pores
- ✓ Brightens dull, uneven skin
- ✓ Minimises pigmentation and sun damage
- ✓ Stimulates collagen for long-term rejuvenation



PLATINUM 3 IN 1 INCLUDES

Picosecond Nd: Yag Laser Tattoo Removal System

Laser Tattoo Removal Treatment

For your clients, Tattoo Removal feels similar to a hot pinprick, whilst this is mildly uncomfortable, most people tolerate treatment well.

- A numbing cream can be applied 20-30 minutes before treatment if required.
- During the treatment, there's an instant whitening of the tattoo, but this fades quickly.
- The area is likely to feel tender for a few hours and appear red and swollen for up to 48 hours.
- With some laser tattoo removal, lasers can kill up to 10% of the hair in the treatment area, though no more than that, and most importantly, often we don't see as much as that.
- Some tattoos will form a graze which will heal naturally.
- Initially, the tattoo is unlikely to look any different but will fade over the next 4 to 8 weeks as the body disperses the damaged ink.
- Very rarely the skin may become lighter or darker, or some blistering may form.
- Another rare side effect is a change to the skin texture and although uncommon, there can be light scarring. It's important to avoid picking grazes or blisters that need to heal naturally. Scarring isn't usually due to the treatment but due to the interference of grazes or blistering.



Tattoo Removal Treatment Steps


Following on the consultation and treatment explanation to the client:

- Disinfect the treatment area.
- Select the right parameters according to the color of the ink.
- Put protective eyewear for yourself and the client.
- Keep the laser head 4-6cm away from the target area and start the treatment.
- It is always recommended to start the treatment with lower energy and gradually increase in each session.
- treatment should be repeated every 4-8 weeks.
- The number of treatments depends on the type and color of the ink, generally 6-8 sessions are recommended.
- Apply Aloe Vera / Repair Gel

Parameters

Ink Colour	Energy j / c m ²	Frequency Hz	Wavelength nm
Black / Blue	600+	1-6	1064
Red / Green	600+	1-6	532
Yellow / Brown	600+	1-6	532

Before

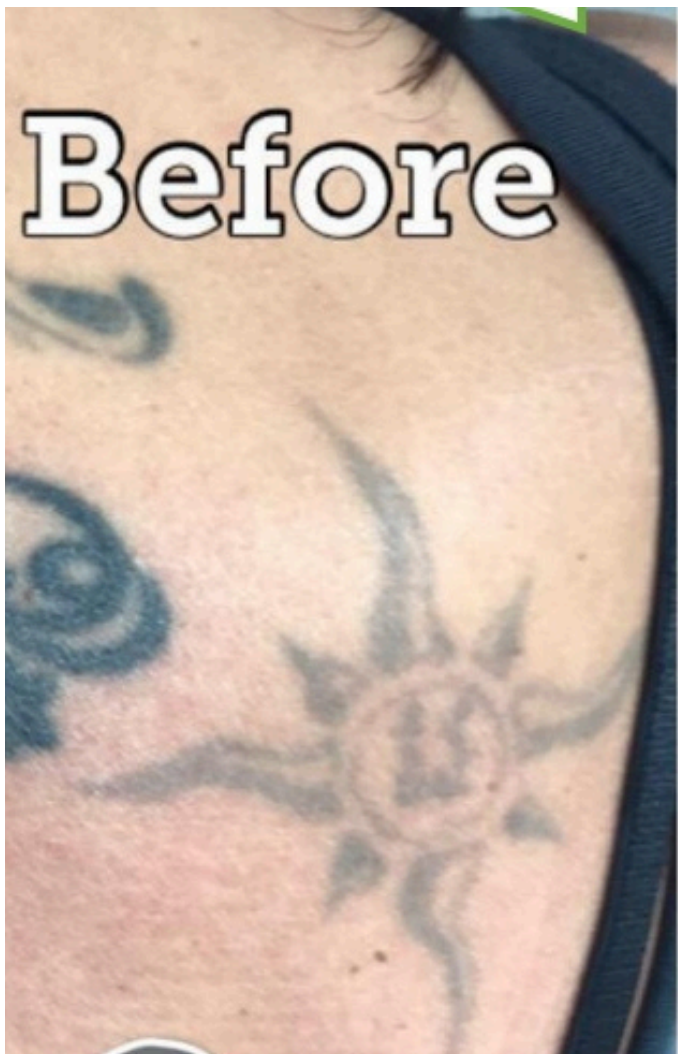

After 1st session




After 3 session



Next level technology
Tattoo removal

PIC•COL

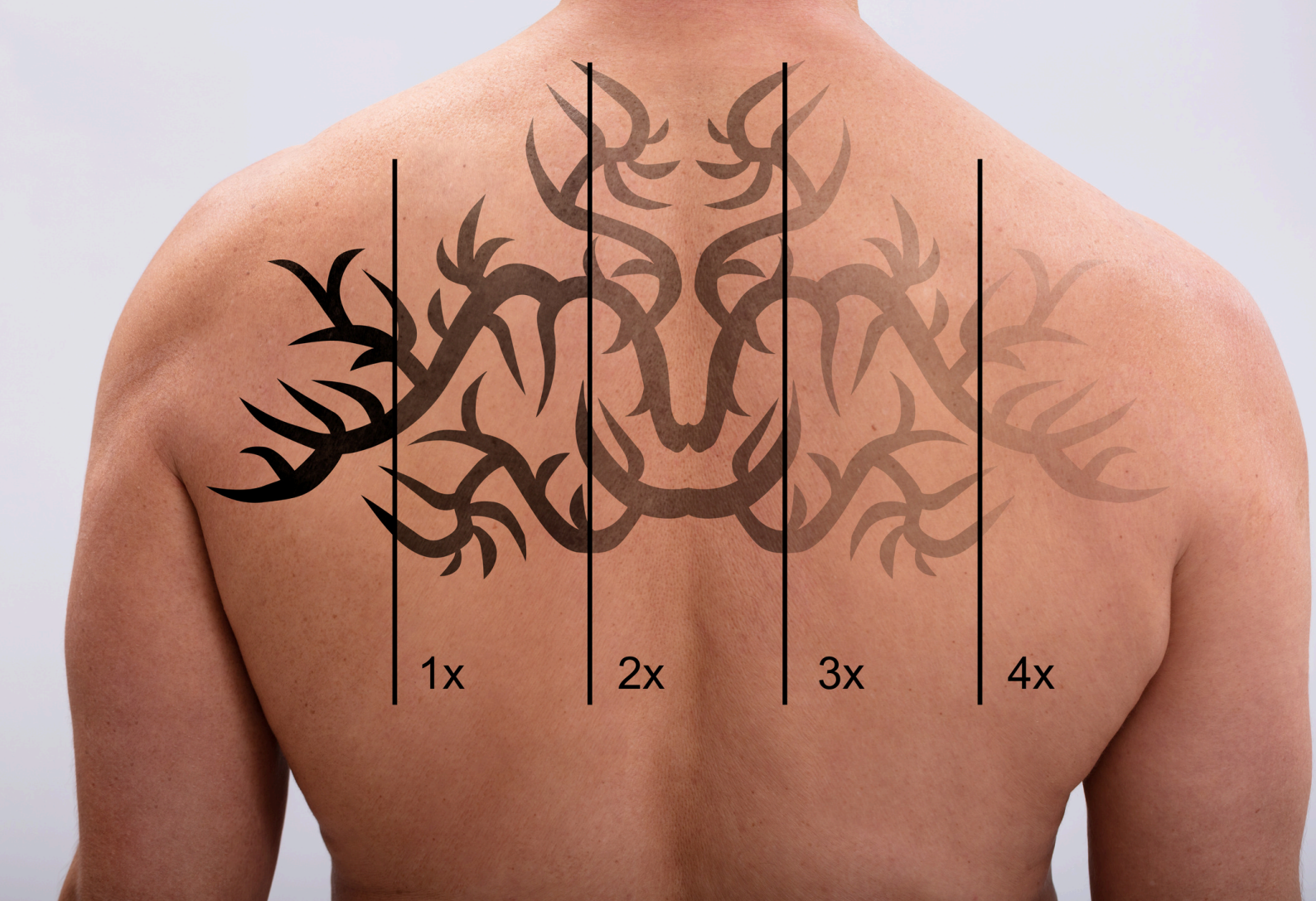

Platinum 3 in 1

L A S E R

**Platinum SkinTech Limited
20-22 Wenlock Road
London
N1 7GU
0203 3627 5009
info@platinumskintechco.uk**